

XNAT ImageViewer: Launch Scenarios

Iteration 2

Author:

Sunil Kumar at Moka Creative, LLC (sunilk@mokacreativellc.com)

Collaborators:

Will Horton at the NRG (hortonw@mir.wustl.edu)

Rick Herrick at the NRG (herrickr@mir.wustl.edu)

Tim Olsen at Deck 5 Consulting, LLC (tim@deck5consulting.com)

ImageViewer: Multi-Subject UI

Multi-Subject Viewer (minimized zippys shown for effect).




- Projects
- + Recent
- + Favorite
- My projects
- Other projects
- + Stored Searches
- + Data

Proj: proj1 +

Subj: subj1 +


Expt: subj1_MR1 +

Scans +

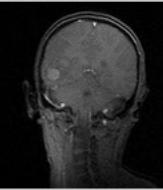

IMA_IGT_GLIOMA_00027_Scan_1_fgre
Frmt: DICOM
Type: LOCALIZER
Expt: XNATImageViewerDemo

Slicer +

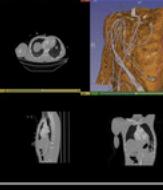

4
Frmt:
Type:
Expt: XNATInr

Expt: subj1_MR2 -


Expt: subj1_MR3 -

Expt: subj1_MR4 -


Subj: subj1a +

Expt: subj1_MR1 +


Scans +

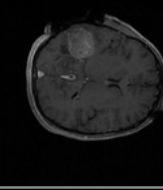

IMA_IC
Frmt: DIC
Type: 3D S
Expt: XNA

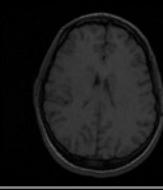


MR-Bra
Frmt: NifTI
Type:
Expt: XNA

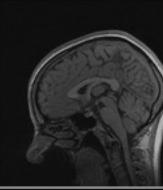

MR-Bra
Frmt: DIC
Type: MPF
Expt: XNA

Multi-Subject Viewer (minimized zippy's shown for effect) with subject spacing.

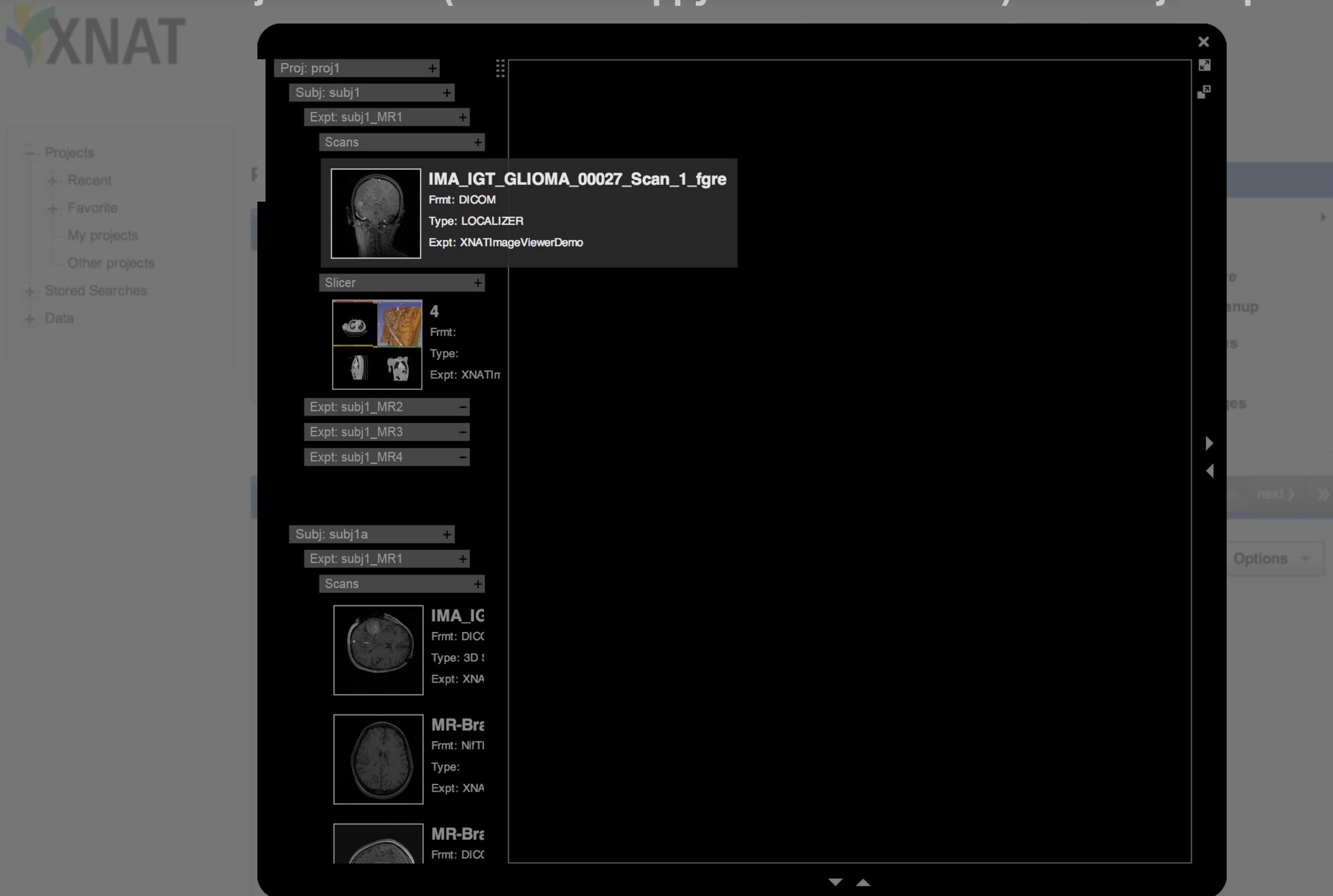

Project Page / Subject Tab

Option A1: “View” links in existing ID column

Logged in as: admin | Auto-logout in: 0 : 14 : 33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Subjects MR Sessions

Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows) Reload Options

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions
subj1 - View	U			4	
subj1a - View	U				1
subj2 - View	U			1	
subj2a - View	U			1	

A red arrow points to the “View” link in the “subj2a” row of the table.

Option A2: “View” links in existing ID column with “Options”

Logged in as: admin | Auto-logout in: 0 : 14 : 33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects **MR Sessions**

Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows)

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions
subj1 - View	U			4	
subj1a - View	U				1
subj2 - View	U			1	
subj2a - View	U			1	

Reload Options

- Spreadsheet
- Email
- Save Search
- Save as New Search
- Show XML
- Edit Columns
- Join to ...
- View Images**

Option B1: “View” column with icons

Logged in as: admin | Auto-logout in: 0 : 14 : 33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects MR Sessions

Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows) Reload Options

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions	View
subj1	U			4		⊕
subj1a	U				1	⊕
subj2	U			1		⊕
subj2a	U			1		⊕

Option B2: “View” column with icons with “Options”

Logged in as: admin | Auto-logout in: 0 : 14 : 33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects MR Sessions

Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows)

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions	View
subj1	U			4		⊕
subj1a	U				1	⊕
subj2	U			1		⊕
subj2a	U			1		⊕

Reload Options

Spreadsheet Email Save Search Save as New Search Show XML Edit Columns Join to ... View Images

Option C: “View” column with check boxes

Logged in as: admin | Auto-logout in: 0 : 14 : 33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

Projects
Recent
Favorite
My projects
Other projects
Stored Searches
Data

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects MR Sessions

Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows) Reload Options

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions	View
subj1	U			4		<input type="checkbox"/>
subj1a	U				1	<input checked="" type="checkbox"/>
subj2	U			1		<input type="checkbox"/>
subj2a	U			1		<input checked="" type="checkbox"/>

[View Selected](#)



Option C: “View” column with check boxes with “Options”

Logged in as: admin | Auto-logout in: 0:14:33 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects **MR Sessions** Add Tab << first < prev 1 next > last >> 20 1 of 1 Pgs (4 Rows) Reload Options

Subject	M/F	Hand	YOB	MR Sessions	CT Sessions	View
subj1	U			4		<input type="checkbox"/>
subj1a	U				1	<input checked="" type="checkbox"/>
subj2	U			1		<input type="checkbox"/>
subj2a	U			1		<input checked="" type="checkbox"/>

View Selected

Options

- Spreadsheet
- Email
- Save Search
- Save as New Search
- Show XML
- Edit Columns
- Join to ...
- View Images

Project Page / MR Sessions Tab

Option A1: “View” links in existing ID column

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Subjects [MR Sessions](#) Add Tab « < prev next > »

<< first < prev 1 next > last >> 20 1 of 1 Pgs (6 Rows) Reload Options

MR ID	Date	Subject	Age	Scanner	Scans
subj1_MR1 - View	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)
subj1_MR2 - View	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)
subj1_MR3 - View	2010-02-14	subj1			Analyze(1)
subj1_MR4 - View	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)
subj2_MR1 - View		subj2			SAG MPRAGE 8 CHANNEL(1)
subj2a_MR1 - View		subj2a			Unknown(1)

A red arrow points to the “- View” link in the Scans column of the second row of the table.

Option A2:

“View” links in existing ID column with “Options”

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects [MR Sessions](#)

Add Tab « < prev next > »

MR ID	Date	Subject	Age	Scanner	Scans
subj1_MR1	- View	2006-01-31	subj1		Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)
subj1_MR2	- View	2006-07-03	subj1		SAG MPRAGE 8 CHANNEL(1)
subj1_MR3	- View	2010-02-14	subj1		Analyze(1)
subj1_MR4	- View	2007-09-04	subj1		SAG MPRAGE 8 CHANNEL(1)
subj2_MR1	- View		subj2		SAG MPRAGE 8 CHANNEL(1)
subj2a_MR1	- View		subj2a		Unknown(1)

Options

- Spreadsheet
- Email
- Save Search
- Save as New Search
- Show XML
- Edit Columns
- Join to ...
- View Images**

Option B1: “View” column with icons

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

[Search](#) [Advanced](#)

[Home](#) [New](#) [Upload](#) [Administer](#) [Tools](#) [Help](#)

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Details [Access](#) [Manage](#) [Pipelines](#)

ID: 1
Description: proj1
Keywords: proj1

[Edit Details](#) [Delete](#) [Manage Custom Variables](#)

Subjects [MR Sessions](#)

Add Tab [«](#) [< prev](#) [next >](#) [»](#)

<< first < prev 1 next > last >> 20 1 of 1 Pgs (6 Rows) [Reload](#) [Options](#)

MR ID	Date	Subject	Age	Scanner	Scans	View
subj1_MR1	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)	
subj1_MR2	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)	
subj1_MR3	2010-02-14	subj1			Analyze(1)	
subj1_MR4	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)	
subj2_MR1		subj2			SAG MPRAGE 8 CHANNEL(1)	
subj2a_MR1		subj2a			Unknown(1)	

[XNATImageViewer: Launch Scenarios Iteration 2] Moka Creative, LLC and the NRG

Option B2:

“View” column with icons with “Options”

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects [MR Sessions](#)

Add Tab « < prev next > »

<< first < prev 1 next > last >> 20 1 of 1 Pgs (6 Rows)

MR ID	Date	Subject	Age	Scanner	Scans	View
subj1_MR1	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)	View
subj1_MR2	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)	View
subj1_MR3	2010-02-14	subj1			Analyze(1)	View
subj1_MR4	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)	View
subj2_MR1		subj2			SAG MPRAGE 8 CHANNEL(1)	View
subj2a_MR1		subj2a			Unknown(1)	View

Options

- Spreadsheet
- Email
- Save Search
- Save as New Search
- Show XML
- Edit Columns
- Join to ...
- View Images

Option C: “View” column with check boxes

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

XNAT

Home New ▾ Upload ▾ Administer ▾ Tools ▾ Help ▾

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage Custom Variables

Subjects MR Sessions

Add Tab « < prev next > »

<< first < prev 1 next > last >> 20 1 of 1 Pgs (6 Rows) Reload Options

MR ID	Date	Subject	Age	Scanner	Scans	View
subj1_MR1	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)	<input type="checkbox"/>
subj1_MR2	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)	<input checked="" type="checkbox"/>
subj1_MR3	2010-02-14	subj1			Analyze(1)	<input type="checkbox"/>
subj1_MR4	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)	<input checked="" type="checkbox"/>
subj2_MR1		subj2			SAG MPRAGE 8 CHANNEL(1)	<input checked="" type="checkbox"/>
subj2a_MR1		subj2a			Unknown(1)	<input type="checkbox"/>

→ View Selected

Option C: “View” column with check boxes with “Options”

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:06 - [renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

proj1

Actions

- Add
- Upload Images
- View Prearchive
- Scan Type Cleanup
- Add to Favorites
- Download XML
- Download Images
- Manage Files

Subjects [MR Sessions](#)

Add Tab « < prev next > »

<< first < prev 1 next > last >> 20 1 of 1 Pgs (6 Rows)

MR ID	Date	Subject	Age	Scanner	Scans	View
subj1_MR1	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)	<input type="checkbox"/>
subj1_MR2	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)	<input checked="" type="checkbox"/>
subj1_MR3	2010-02-14	subj1			Analyze(1)	<input type="checkbox"/>
subj1_MR4	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)	<input checked="" type="checkbox"/>
subj2_MR1		subj2			SAG MPRAGE 8 CHANNEL(1)	<input type="checkbox"/>
subj2a_MR1		subj2a			Unknown(1)	<input type="checkbox"/>

Options

- Spreadsheet
- Email
- Save Search
- Save as New Search
- Show XML
- Edit Columns
- Join to ...
- View Images

View Selected

Subject Page

Logged in as: admin | ? Auto-logout in: 0:14:51 - [renew](#) | [Logout](#)

Home New Upload Administer Tools Help

 Search Advanced

- Projects
- + Recent
- + Favorite
- ... My projects
- ... Other projects
- + Stored Searches
- + Data

[PROJECT: proj1 > subj1](#)

Subject Details: subj1

Details Projects

Accession # localhost_S00002
Date Added 2013-12-28 15:13:46.622 (admin)
Birth year --
Gender
Handedness

Actions

Edit
View XML
Add Experiment
Upload Images
Download XML
Email
Manage Files
Delete

Experiments

Date	Experiment	Project	Label
2006-01-31	MR Session - View	proj1	subj1_MR1
2006-07-03	MR Session - View	proj1	subj1_MR2
2007-09-04	MR Session - View	proj1	subj1_MR4
2010-02-14	MR Session - View	proj1	subj1_MR3



Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:51 - [renew](#) | [Logout](#)

XNAT version 1.6.2.1

Home New Upload Administer Tools Help

PROJECT: proj1 > subj1

Subject Details: subj1

Details Projects

Accession # localhost_S00002
Date Added 2013-12-28 15:13:46.622 (admin)
Birth year --
Gender
Handedness

Actions

Edit
View XML
Add Experiment
Upload Images
Download XML
Email
Manage Files
Delete

Experiments

Date	Experiment	Project	Label	View
2006-01-31	MR Session	proj1	subj1_MR1	
2006-07-03	MR Session	proj1	subj1_MR2	
2007-09-04	MR Session	proj1	subj1_MR4	
2010-02-14	MR Session	proj1	subj1_MR3	

Logged in as: [admin](#) | [?](#) Auto-logout in: 0:14:51 - [renew](#) | [Logout](#)

XNAT
XNAT version 1.6.2.1

Home New Upload Administer Tools Help

PROJECT: proj1 > subj1

Subject Details: subj1

Details Projects

Accession # localhost_S00002
Date Added 2013-12-28 15:13:46.622 (admin)
Birth year --
Gender
Handedness

Actions

Edit
View XML
Add Experiment
Upload Images
Download XML
Email
Manage Files
Delete

Experiments

Date	Experiment	Project	Label	View
2006-01-31	MR Session	proj1	subj1_MR1	<input type="checkbox"/>
2006-07-03	MR Session	proj1	subj1_MR2	<input checked="" type="checkbox"/>
2007-09-04	MR Session	proj1	subj1_MR4	<input type="checkbox"/>
2010-02-14	MR Session	proj1	subj1_MR3	<input checked="" type="checkbox"/>

[View Selected](#)



[XNATImageViewer: Launch Scenarios Iteration 2]

Moka Creative, LLC and the NRG

Experiment Page

Actions option and “View” links in inner column

XNAT

Logged in as: admin | Auto-Logout: 0:14:47 - [view] | Layout
 Search Advanced

Home New Upload Administer Tools Help

PROJECT: proj1 > SUBJECT: subj1 > subj1_MR1

MR Session: subj1_MR1

Active Processes

- AutoRun: Queued Start Time: 2013-12-31 08:13:03.916 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-31 08:05:28.25 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-31 07:58:19.999 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-29 07:41:25.295 [Mark as Failed]

Details	Projects	Actions
Accession #: localhost_E00007 Date Added: 12/29/2013 07:41:24 (admin) Date: 01/31/2006 Time: 09:18:45 Scanner: GE MEDICAL SYSTEMS GENESIS_SIGNA Acquisition Site: Brain Mapping Center	Subject: subj1 Gender: Handedness: Age: -	Edit View ► Upload ► Download ► Email View in XImgView Manage Files Delete

Scans

Scan	Type	View	Series Desc	Usability	Files	Note
+ IMA	Unknown	View	SAG MPRAGE 8 CHANNEL	usable	7.8 MB in 63 files	
+ 6	SAG MPRAGE 8 CHANNEL	View	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files	
+ 7	SAG MPRAGE 8 CHANNEL	View	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files	
+ 8	Analyze	View		unsupported	24.0 MB in 2 files	

Unsupported formats do not have the option

Total: 72.5 MB in 385 files

Actions option and “View” links in new column



Logged in as: admin | Auto-logout in 0:41:47 - [Logout](#)

[Search](#) [Advanced](#)

Home New Upload Administer Tools Help

PROJECT: proj1 > SUBJECT: subj1 > subj1_MR1

MR Session: subj1_MR1

Active Processes	
AutoRun: Queued Start Time: 2013-12-31 08:13:03.916	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-31 08:05:28.25	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-31 07:58:19.999	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-29 07:41:25.295	[Mark as Failed]

Details **Projects**

Accession #: localhost_E00007 **Subject:** subj1
Date Added: 12/29/2013 07:41:24 (admin) **Gender:**
Date: 01/31/2006 **Handedness:**
Time: 09:18:45 **Age:** —
Scanner: GE MEDICAL SYSTEMS GENESIS_SIGNA
Acquisition Site: Brain Mapping Center

Actions

- Edit
- View
- Upload
- Download
- Email
- View in XImgView
- Manage Files
- Delete

A red arrow points from the "View" link in the Actions column to the "View" link in the "Scans" table row for the SAG MPRAGE 8 CHANNEL series.

Scans

Scan	Type	Series Desc	Usability	Files	Note	View
+ IMA	Unknown		usable	7.8 MB in 63 files		View
+ 6	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		View
+ 7	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		View
+ 8	Analyze		usable	24.0 MB in 2 files		

Unsupported formats do not have the option

Actions option and “View” icons in new column



Logged in as: admin | Auto-Logout in: 0:14:47 - Logout | Search | Advanced

Home New Upload Administer Tools Help

PROJECT: proj1 > SUBJECT: subj1 > subj1_MR1

MR Session: subj1_MR1

Active Processes

- AutoRun: Queued Start Time: 2013-12-31 08:13:03.916 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-31 08:05:28.25 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-31 07:58:19.999 [Mark as Failed]
- AutoRun: Queued Start Time: 2013-12-29 07:41:25.295 [Mark as Failed]

Details	Projects	Actions
Accession #: localhost_E00007 Date Added: 12/29/2013 07:41:24 (admin) Date: 01/31/2006 Time: 09:18:45 Scanner: GE MEDICAL SYSTEMS GENESIS_SIGNA Acquisition Site: Brain Mapping Center	Subject: subj1 Gender: Handedness: Age: -	Edit View ➔ Upload ➔ Download ➔ Email View in XImgView Manage Files Delete

Scans

Scan	Type	Series Desc	Usability	Files	Note	View
+ IMA	Unknown	SAG MPRAGE 8 CHANNEL	usable	7.8 MB in 63 files		🕒
+ 6	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		🕒
+ 7	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		🕒
+ 8	Analyze		usable	24.0 MB in 2 files		

Actions option and “View” check boxes in new column

XNAT

Home New Upload Administer Tools Help

PROJECT: proj1 > SUBJECT: subj1 > subj1_MR1

MR Session: subj1_MR1

Active Processes	
AutoRun: Queued Start Time: 2013-12-31 08:13:03.916	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-31 08:05:28.25	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-31 07:58:19.999	[Mark as Failed]
AutoRun: Queued Start Time: 2013-12-29 07:41:25.295	[Mark as Failed]

Details Projects

Accession #: localhost_E00007 Subject: subj1
 Date Added: 12/29/2013 07:41:24 (admin) Gender:
 Date: 01/31/2006 Handedness:
 Time: 09:18:45 Age:
 Scanner: GE MEDICAL SYSTEMS GENESIS_SIGNA
 Acquisition Site: Brain Mapping Center

Actions

- Edit
- View
- Upload
- Download
- Email
- View in XImgView
- Manage Files
- Delete



Scans

Scan	Type	Series Desc	Usability	Files	Note	View
+ IMA	Unknown		usable	7.8 MB in 63 files		<input checked="" type="checkbox"/>
+ 6	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		<input checked="" type="checkbox"/>
+ 7	SAG MPRAGE 8 CHANNEL	SAG MPRAGE 8 CHANNEL	usable	20.4 MB in 160 files		<input checked="" type="checkbox"/>
+ 8	Analyze		usable	24.0 MB in 2 files		



Unsupported formats do not have the option

“Load into open viewer?” Popup

Logged in as: admin | [Auto-logout in: 0:14:06 - renew](#) | [Logout](#)

XNAT

Home New Upload Administer Tools Help

Projects
Recent
Favorite
My projects
Other projects
Stored Searches
Data

proj1

Details Access Manage Pipelines

ID: 1
Description: proj1
Keywords: proj1

Edit Details Delete Manage C View

Add
Upload Images
View Prearchive
Scan Type Cleanup

Add to Favorites
Download XML
Download Images
Manage Files

View

Load into open viewer?

Subjects MR Sessions


<< first < prev 1 next> last>> 20

Yes No Reload Options


MR ID	Date	Subject	Age	Scanner	Scans	View
subj1_MR1	2006-01-31	subj1			Analyze(1), SAG MPRAGE 8 CHANNEL(2), Unknown(1)	
subj1_MR2	2006-07-03	subj1			SAG MPRAGE 8 CHANNEL(1)	
subj1_MR3	2010-02-14	subj1			Analyze(1)	
subj1_MR4	2007-09-04	subj1			SAG MPRAGE 8 CHANNEL(1)	
subj2_MR1		subj2			SAG MPRAGE 8 CHANNEL(1)	
subj2a_MR1		subj2a			Unknown(1)	


Appendix: Sketches

Project Page / Subject Tab

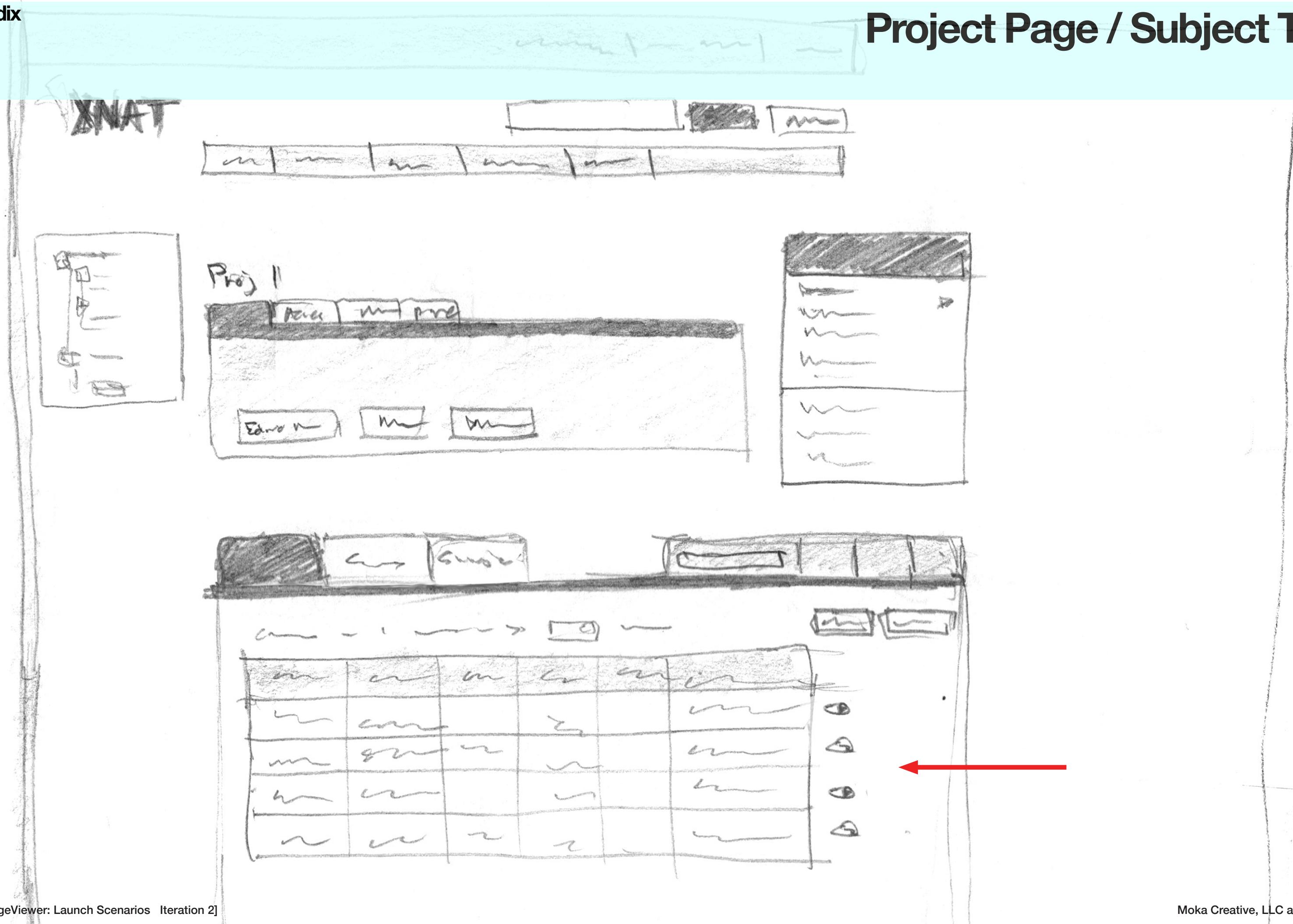

Project Page / MR Sessions Tab

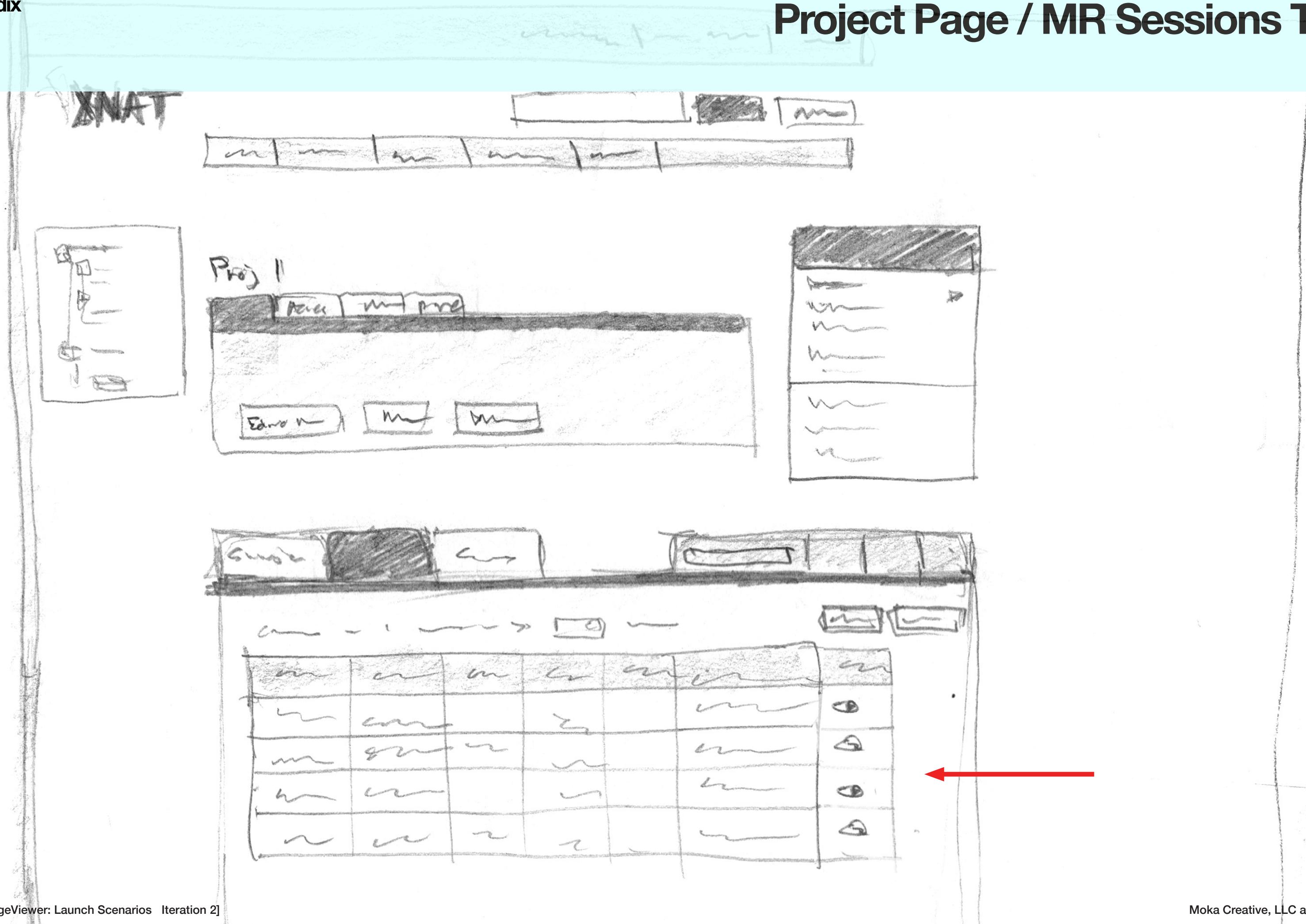


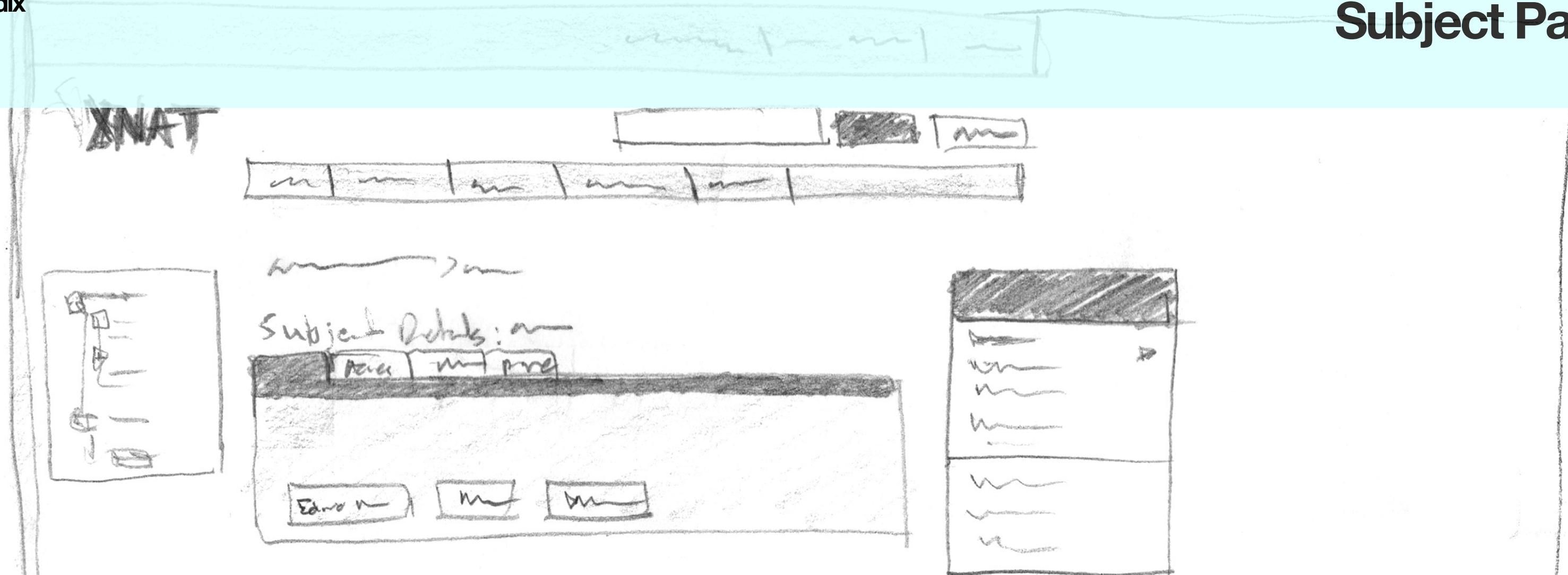

Experiments

Date	Entered	Project	Last
w	w	w	w
w	w	w	w
w	w	w	w
w	w	w	w



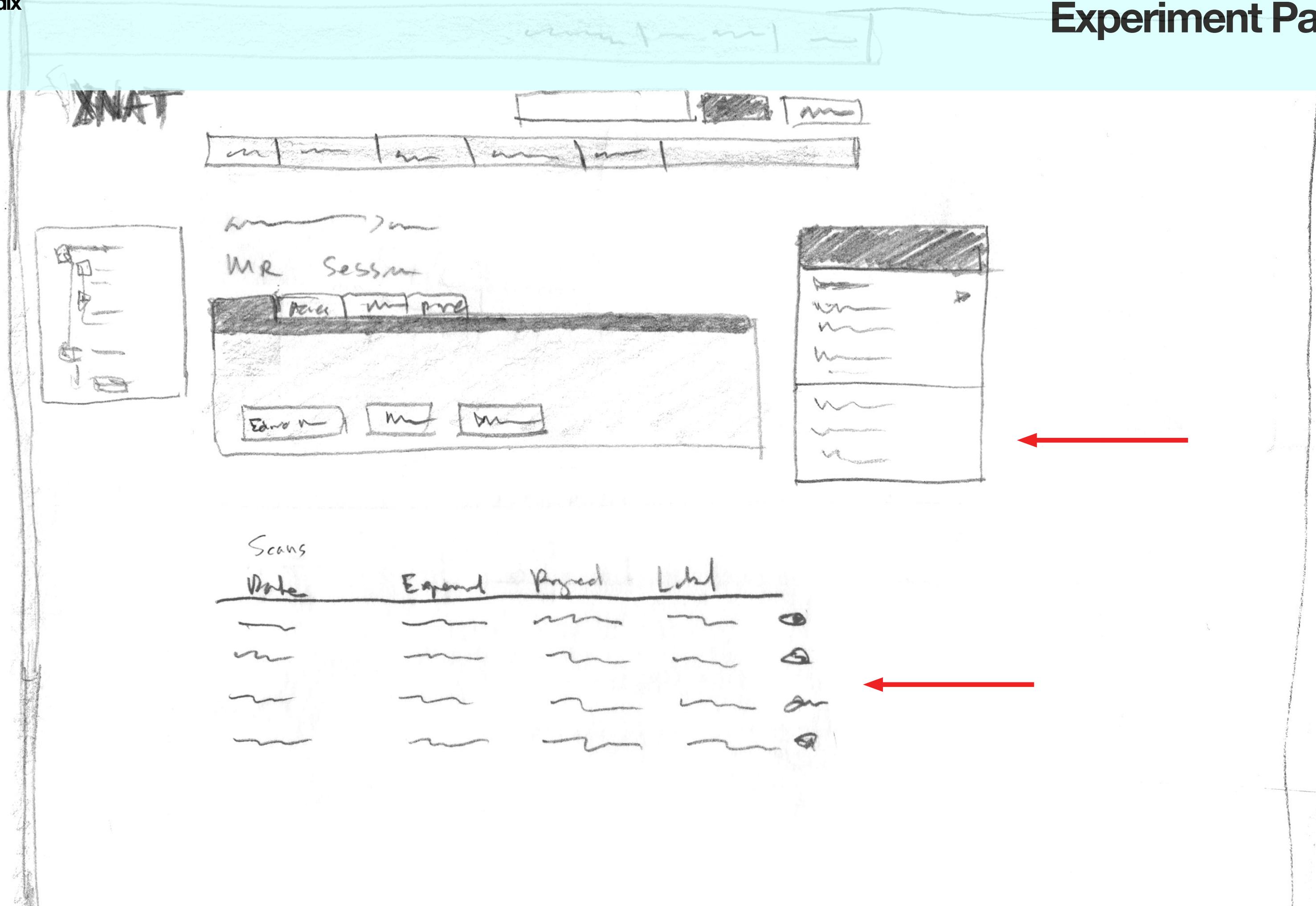

“Load into open” popup

